

# Body Shape RF

## Fractional Micro-Needles & Quadro Polar Technologies

**BODY SHAPE RF** is an innovative product, designed to achieve:

- Skin Rejuvenation
- Body Contouring
- Skin Lifting
- Cellulite Improvement

### INNOVATIVE TECHNOLOGY:

- **Temperature Sensors** integrated in the Applicators for continuous skin-temperature monitoring, assuring effective and problem free treatments
- **Fractional Micro-Needles** for effective and deep Skin Rejuvenation and Facial Wrinkle Removal.
- **Quadro & Bi-Polar Applicators:** Large working surface, with exact layer and deep tissue targeting
- **High Power** (up to 250 W<sub>p</sub>) enabling a powerful treatment with the best possible results



# Body Shape RF

The Safe Method for Body Contouring

## Bi-Polar



Temp. Controlled  
Facial Applicator



• Deep Energy •  
Penetration

## Quadro-Polar

Temp. Controlled  
Body Applicator

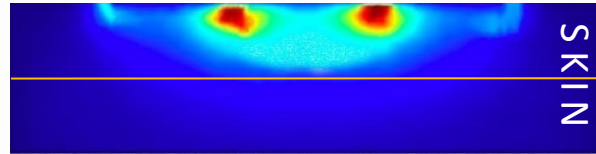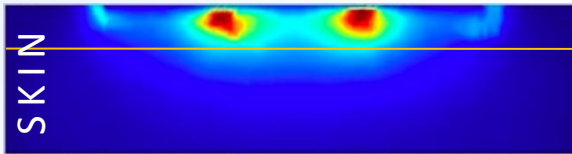


## Physiological Process

Twisted & unresponsive  
Collagen Fibers


Visible Wrinkle

Sorted & supportive  
Collagen Fibers

Flattened Wrinkle

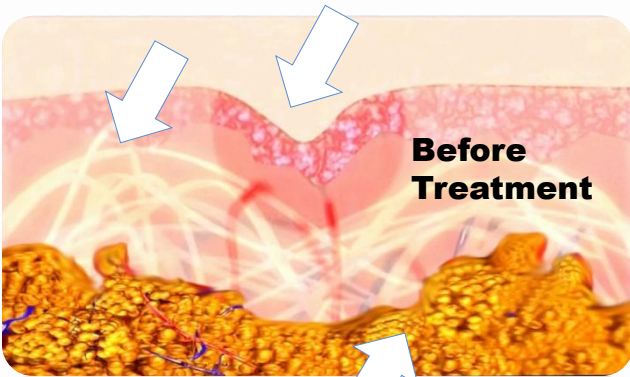

Excessive layer of  
Enlarged Fat Cells

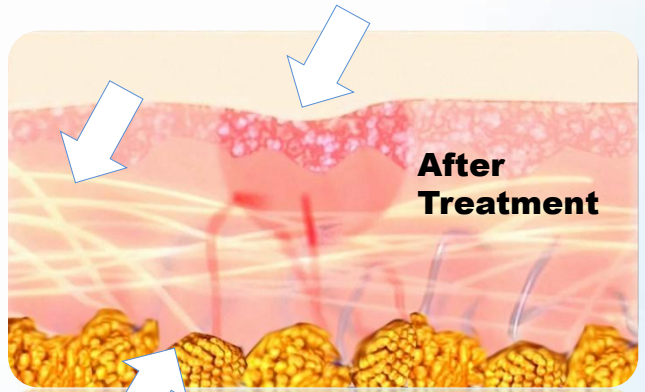


Reduced layer of  
Normal Fat Cells

# Body Shape RF


## Fractional Micro-Needles & Quadro Polar Technologies

### *Fractional RF* with Disposable **Micro-Needles** 6x6 Heads


A built-in Fractional Handpiece with a selection of 3 types of Micro Needles Heads (0.5, 1.0 and 1.5 mm), is the perfect tool for Facial Skin rejuvenation and Fine Wrinkle Removal. The tiny “Holes” are triggering the skin to heal the “damage” by replacing cells and generating Collagen Fibers.

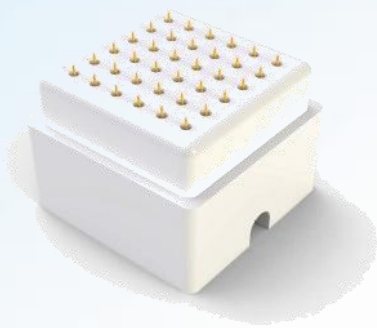

During RF Emission

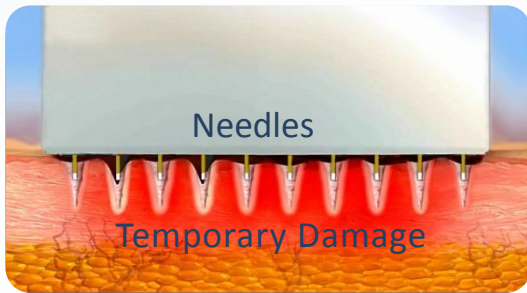

After Removing Needles

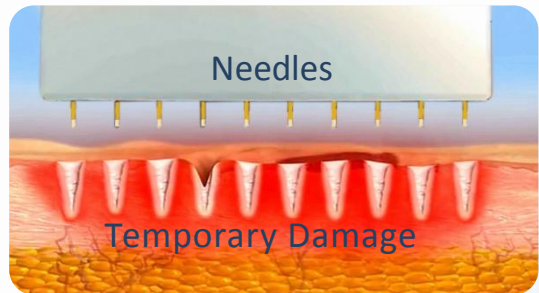

Before



After 2 Treatments



# Body Shape RF

Fractional Micro-Needles & Quadro Polar Technologies



## 3 Simple Steps Operation:

- a. Set the Time Duration
- b. Set the recommended Energy
- c. Set the targeted Skin Temperature
- D. On-Line Temperature Detector**

Before



After



Before



After



MAXIMUM OUTPUT	250 W <sub>P</sub> (100W RMS)
FREQUENCY	6.8 MHZ
APPLICATORS	BODY QUADRO POLAR FACIAL BI-POLAR
TEMPERATURE SENSOR	INFRA-RED HI-SPEED INTEGRATED TEMPERATURE SENSOR IN EACH APPLICATOR
EMERGENCY HALTING	PATIENT'S EMERGENCY HALTING KNOB
MAINS	220-230VAC 50/60HZ 100-120VAC 50/60HZ
MAX MAINS POWER	600 VA